

AD A 044844

NSWC/WOL/TR 77-71


12

ANALYSIS OF MERCURIC 5-NITROTETRAZOLE

BY DONALD J. GLOVER

RESEARCH AND TECHNOLOGY DEPARTMENT

19 JULY 1977

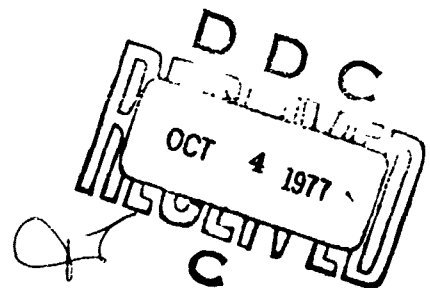

AD No. _____
DDC FILE COPY



NAVAL SURFACE WEAPONS CENTER

Dahlgren, Virginia 22448 • Silver Spring, Maryland 20910

REPRODUCTION STATEMENT A
Approved for public release;
Distribution Unlimited

UNCLASSIFIED

SECURITY CLASSIFICATION OF THIS PAGE (When Data Entered)

| REPORT DOCUMENTATION PAGE | | READ INSTRUCTIONS BEFORE COMPLETING FORM |
|---|---|---|
| 1. REPORT NUMBER NSWC/WOL/TR-77-71 | 2. GOVT ACCESSION NO. | 3. RECIPIENT'S CATALOG NUMBER |
| 4. TITLE (and Subtitle) ANALYSIS OF MERCURIC 5-NITROTETRAZOLE | | 5. TYPE OF REPORT & PERIOD COVERED Final Rept. |
| 6. AUTHOR DONALD J. GLOVER | | 7. PERFORMING ORG. REPORT NUMBER |
| 8. PERFORMING ORGANIZATION NAME AND ADDRESS Naval Surface Weapons Center White Oak Laboratory White Oak, Silver Spring, Maryland 20910 | | 9. PROGRAM ELEMENT, PROJECT, TASK AREA & WORK UNIT NUMBERS OPN; 0; 0; WRL1HA; |
| 10. CONTROLLING OFFICE NAME AND ADDRESS | 11. REPORT DATE 19 July 1977 | 12. NUMBER OF PAGES 12 |
| 13. MONITORING AGENCY NAME & ADDRESS (if different from Controlling Office) 12-16f. | 14. SECURITY CLASS (of this report) UNCLASSIFIED | |
| 15a. DECLASSIFICATION/DOWNGRADING SCHEDULE | | |
| 15. DISTRIBUTION STATEMENT (of this Report) Approved for public release; distribution unlimited | | |
| 17. DISTRIBUTION STATEMENT (of the abstract entered in Block 20, if different from Report) | | |
| 18. SUPPLEMENTARY NOTES | | |
| 19. KEY WORDS (Continue on reverse side if necessary and identify by block number) Mercury (II) UV absorption 5-nitrotetrazole dichromate oxidation primary explosives thiocyanate titration analyses atomic absorption | | |
| 20. ABSTRACT (Continue on reverse side if necessary and identify by block number) Mercury (II) in mercuric 5-nitrotetrazole may be determined by atomic absorption using the unique-to-mercury absorption line at 2536.5 angstroms, with a precision of 1-3%. For a precision of 0.1 to 0.6%, a titration with thiocyanate ion is recommended. The 5-nitrotetrazole anion may be determined by measuring the absorbance due to the anion at 257 nm in the ultraviolet, with a precision of + 1%. The anion may also be determined by oxidation with dichromate ion in sulfuric acid solution, with a precision | | |

DD FORM 1 JAN 73 1473

EDITION OF 1 NOV 65 IS OBSOLETE
S/N 0100-107-014-6001

UNCLASSIFIED

SECURITY CLASSIFICATION OF THIS PAGE (When Data Entered)

391 596

JP

UNCLASSIFIED

SECURITY CLASSIFICATION OF THIS PAGE (When Data Entered)

→ of $\pm 0.5\%$. ←

$\pm 0.5\%$

UNCLASSIFIED

SECURITY CLASSIFICATION OF THIS PAGE (When Data Entered)

NSWC/WOL/TR 77-71

19 July 1977

ANALYSIS OF MERCURIC 5-NITROTETRAZOLE

This report describes analytical procedures for the determination of the mercury (II) cation, as well as the 5-nitrotetrazole anion, in mercuric 5-nitrotetrazole, in particular for use in a specification. The work was carried out under Task No. WR11HA(OPN).

Julius W. Enig

JULIUS W. ENIG
By direction

| | |
|----------------------------------|-------------------------------------|
| COPIES FOR | |
| White Section | <input checked="" type="checkbox"/> |
| Blue Section | <input type="checkbox"/> |
| BY | |
| DISTRIBUTION/AVAILABILITY COPIES | |
| BY CIAL | |
| A | |

TABLE OF CONTENTS

| | Page |
|------------------------------------|------|
| INTRODUCTION | 3 |
| EXPERIMENTAL AND RESULTS | 4 |
| DISCUSSION | 9 |
| SUMMARY AND CONCLUSIONS | 11 |




| Table | TABLES Title | Page |
|-------|---|------|
| I | Mercuric 5-Nitrotetrazole, Anion Determination | 5 |
| II | Mercuric 5-Nitrotetrazole, Determination of Mercury and Copper | 7 |

ANALYSIS OF MERCURIC 5-NITROTETRAZOLE

INTRODUCTION

Mercuric 5-nitrotetrazole ($\text{Hg}(5\text{-NT})_2$) is being studied at the NSWC, WOL as a possible replacement for lead azide (1). Unlike lead azide, $\text{Hg}(5\text{-NT})_2$ is stable to moisture and is compatible with detonator construction materials, such as copper or brass.


Twelve batches of $\text{Hg}(5\text{-NT})_2$ ranging from 65 g to 600 g were recently prepared (2) by the procedure of Gilligan and Kamlet (1), according to the following series of reactions:

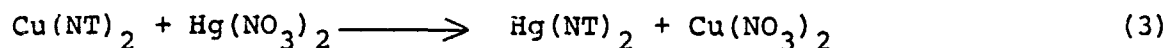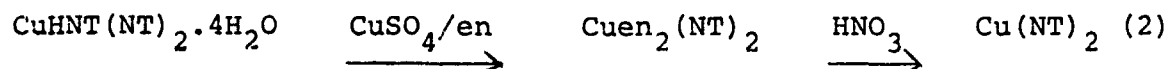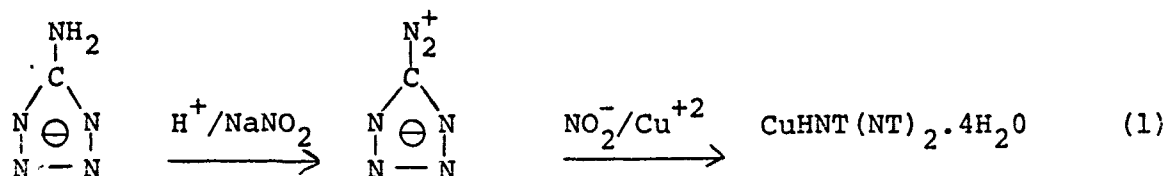

where $\text{NT} = \text{NO}_2$, the nitrotetrazole anion, and $\text{en} = \text{H}_2\text{NCH}_2\text{CH}_2\text{NH}_2$,



ethylenediamine. In addition to these desired reactions, the initial diazo compound (equation (1)) can also react to form the hydroxy-tetrazole (1):

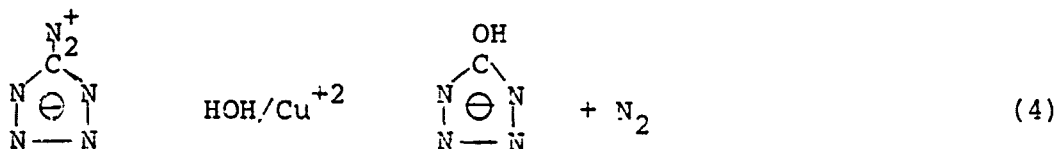

1. W. H. Gilligan and M. J. Kamlet, "Synthesis of Mercuric 5-nitrotetrazole", NSWC/WOL/TR 76-146, 9 December 1976.

2. R. E. Farncomb, M. Chang, and F. J. Pisacane, "Process Scale-up for Mercuric 5-nitrotetrazole", personal communication.

It is desirable to have analytical methods for both Hg(II) and the anion, NT, in the final product, as well as for any impurities.

A method is already described (3) for the anion, based on the maximum absorption in the ultraviolet at 257nm. This method is rapid and works well. Another method for the anion, oxidation with dichromate in sulfuric acid, is presented here.

Of the numerous methods in the literature for Hg(II), two were selected for evaluation. A method for determining Hg(II) by atomic absorption (AA) would confirm that this is the element being measured. For this reason, the AA method presented in this report is recommended for specification purposes, even though the method gives a precision of only 1 to 3%. Greater precision may be obtained by the thiocyanate procedure described here.

EXPERIMENTAL AND RESULTS

CAUTION: ANY SALT OF 5-NITROTETRAZOLE SHOULD BE CONSIDERED HAZARDOUS WHEN HANDLED. THE MERCURY (II) SALT HAS THE SAME IMPACT SENSITIVITY, 5 CM., EITHER WET OR DRY (reference 2).

Determination of 5-nitrotetrazole anion

UV Method

Accurately weigh approximately 2.5 mg of Hg (5-NT)₂ into a 25 ml volumetric flask. Dissolve in 2.5 ml of 20% ammonium acetate solution, dilute to volume with water, and mix. This is Solution A.

Pipet a two ml aliquot of Solution A into a 10 ml volumetric flask, dilute to volume with water, mix, and read the absorbance at 257nm. Frequently, acetone is used in cleaning volumetric glassware. Even a trace of acetone will interfere at 257 nm and must not be used with any glassware used in this determination.

The ammonium acetate solution is necessary for dissolution of the Hg (5-NT)₂ in a reasonable time. Hg (5-NT)₂ is soluble in water alone, but slowly, 89.7 mg/100ml dissolving in 3 days at 30°, with 95% of this dissolving in 7 hours.

Using the above procedure with a specially purified sample of Hg (5-NT)₂, the molar absorbancy was found to be 5353 ± 28 (average and average deviation for four samples). NMR was used to establish

(3) N. J. Blay, D. G. Davies, D. C. Mullenger, and R. J. Rapley, "Silver Nitrotetrazole: Its Stability and Compatibility with Other Materials", ERDE Technical Report No. 163, January 1974.

the absence of water and any material having hydrogen atoms. A sample of sodium 5-NT was recrystallized from acetone and air dried. If dried in a desiccator, this salt lost its water of hydration (NMR), but if left loosely stoppered on the lab bench, it was shown by NMR to be the di-hydrate. Four samples of this di-hydrate gave 5329 ± 10 for the molar absorptancy of the NT anion at 257nm.

Fourteen different samples of Hg (5-NT)₂ were analyzed by the UV method and the results are shown in Table I.

Dichromate method (4)

Accurately weigh a 35-44 mg sample of Hg (5-NT)₂ and transfer it to a 250 ml standard tapered Erlenmeyer flask, using 20 ml of water. Pipet 10 ml of 0.2500 N K₂Cr₂O₇ into the flask, then add 30 ml of 96% sulfuric acid, which contains 10 g of Ag₂SO₄/liter. Mix thoroughly. Connect the flask to a reflux condenser and reflux at least 30 minutes. Carefully rinse the condenser with distilled water, and cool the flask to room temperature. Dilute the solution to about 125 ml, again cool, add 2 drops of ferroin indicator and titrate with standardized 0.1 N ferrous ammonium sulfate solution. The end point is very sharp, 0.02 ml, going from green to light brown.

All glassware must be cleaned by using the above procedure without the sample. The distilled water must not be stored in plastic wash bottles. A blank must be run using everything except the sample. The milliequivalents of dichromate consumed are calculated by subtracting the sample titre from the blank titre, and multiplying the difference by the normality of the Fe⁺². When the milliequivalents thus found are divided by the millimoles of anion (= 2 x millimoles of Hg (5-NT)₂), a value of 6 is found, which corresponds to the expected electron change (see Discussion).

The Ag₂SO₄ is used as a catalyst (4), and if it is omitted, the reaction is only 80-90% complete in 30 minutes. With the catalyst, the reaction is complete in 15 minutes.

The results using this procedure are given in Table I, for both Hg (5-NT)₂ and Na 5-NT.2HOH.

Determination of Mercury (II) cation

AA Method

Prepare a standard solution of 1000 ppm of Hg(II) by dissolving 1.080 g of mercuric oxide in 5 ml of 6 N nitric acid in a one liter volumetric flask, and diluting to volume with water. Do not use

(4) "Standard Methods for the Examination of Water and Wastewater", edited by H. P. Orland, 12th ed. 1965, p. 510, "Chemical Oxygen Demand".

TABLE I
MERCURIC 5-NITROTETRAZOLE ANION DETERMINATION

| Sample (a) | UV Absorption | | | | Dichromate Oxidation | | | | | |
|--------------|------------------|------------------|--------------|-----------------------------|----------------------|------------------|---------------|--------------|-----------------------------|--------------|
| | No. of Detns. | Percent Found | Ave. Dev. | Percent of calcd. (b) | Ave. Dev. | No. of Detns. | Meq. Mmole | Ave. Dev. | Percent of calcd. (c) | Ave. Dev. |
| 96-8898-80 | 5 | 53.37 | 0.32 | 100.3 | 0.6 | 10 | 5.983 | 0.025 | 99.7 | 0.42 |
| Modified | 6 | 51.83 | 0.21 | 97.4 | 0.4 | 2 | 6.045 | 0.015 | 100.8 | 0.25 |
| WOL(M) | 8 | 52.62 | 0.32 | 98.9 | 0.6 | 2 | 6.034 | 0.030 | 100.6 | 0.50 |
| WOL(P) | 6 | 52.52 | 0.53 | 98.7 | 1.0 | 4 | 6.042 | 0.025 | 100.7 | 0.04 |
| Batch 2 | 7 | 52.47 | 0.42 | 98.6 | 0.8 | 2 | 6.019 | 0.032 | 100.3 | 0.53 |
| " 3 | 8 | 52.31 | 0.42 | 98.3 | 0.8 | 2 | 6.098 | 0.026 | 101.6 | 0.43 |
| " 4 | 7 | 52.84 | 0.58 | 99.3 | 1.1 | 2 | 5.990 | 0.016 | 99.8 | 0.27 |
| " 5 | 6 | 52.52 | 0.32 | 98.7 | 0.6 | 2 | 5.990 | 0.020 | 99.8 | 0.33 |
| " 6 | 9 | 47.46 | 0.28 | 89.2 | 0.6 | 2 | 5.649 | 0.028 | 94.2 | 0.50 |
| " 7 | 6 | 52.47 | 0.52 | 98.6 | 1.0 | 1 | 5.975 | -- | 99.6 | -- |
| " 9 | 6 | 52.73 | 0.47 | 99.1 | 0.9 | 2 | 5.964 | 0.005 | 99.4 | 0.08 |
| " 10 | 5 | 53.10 | 0.32 | 99.8 | 0.6 | 2 | 5.986 | 0.020 | 99.8 | 0.33 |
| " 11 | 5 | 52.89 | 0.32 | 99.4 | 0.6 | 2 | 6.000 | 0.019 | 100.0 | 0.32 |
| " 12 | 5 | 52.78 | 0.53 | 99.2 | 1.0 | 2 | 5.930 | 0.028 | 98.8 | 0.47 |
| Na 5-NT.2HOH | | | | | | 4 | 5.988 | 0.012 | 99.8 | 0.20 |

(a) The first sample is a laboratory prep.; all others are large scale, 65g to 600g batches; WOL(M) and WOL(P) are either Batch 1 or 8 of a series of 12 batches.

(b) Calcd. for Hg (5-NT)₂; anion = 53.21%; Hg = 46.79%.

(c) Based on 6 electron change.

Hg(I) compounds as standards, as they give erroneously high results (5). Dilute aliquots of this standard to give solutions containing 30, 40, and 50 ppm of Hg(II). Read the absorbancy of Solution A above at 2536.5 angstroms and compare with the diluted standards. The machine sensitivity must be adjusted so that the 50 ppm solution reads about three-fourths of the dial. For the Varian Techtron AA-5, this reading is 75 using the times 10 control.

The results using this procedure on Hg (5-NT)₂ are given in Table II.

Thiocyanate Method

Pipet an aliquot of Solution A into a suitable titration vessel (a 25 ml volumetric flask may be used), add 1 ml of 6 N nitric acid, 1 ml of 0.1 N ferric nitrate, mix and titrate with 1×10^{-3} M sodium thiocyanate solution to the light brown end point. Determination of the end point is facilitated by comparison with the color of the starting solution and a solution of the standard Hg(II) that has been titrated. A solution that is 50 ppm in Hg(II) is 2.5×10^{-4} M. The end point occurs when two moles of SCN⁻ per mole of Hg(II) have been added. Consequently, an aliquot of 8 ml will require 4 ml of SCN⁻. A micro buret graduated to 0.01 ml is therefore required.

The results using this procedure are given in Table II.

The end point is easier to detect, if about 50 mg of Hg (5-NT)₂ are titrated using 0.1 N SCN⁻. The sample is dissolved in 2.5 ml of 20% ammonium acetate, and because ferric ion gives a red color with ammonium ion, 3 ml of 6 N HNO₃ are required to discharge this color. The more dilute solution above was used to avoid unnecessary handling of sensitive Hg (5-NT)₂.

Determination of Copper (II) cation

AA Method

Copper (II) is a likely impurity in Hg (5-NT)₂, because it is used during the synthesis. A 1 ppm solution of Cu(II) gives a dial reading of 100 on the times 10 dial setting on the Varian Techtron AA-5, whereas 50 ppm Hg(II) gives a dial reading of 75. Therefore, up to 2% Cu(II) in Hg(II) can be determined using the same solution. On this instrument, the Cu and Hg lamps may be mounted at the same time. After determining the absorption due to Hg(II), one changes the lamp by rotating a disk, adjusts the wavelength, and then determines the absorption due to Cu(II).

The results using this procedure are given in Table II.

(5) D. N. Hingle, G. F. Kirkbright, and T. S. West, *Analyst*, 92, 759 (1967); *Chem. Abstracts*, 68, 111,101g (1968).

TABLE II
MERCURIC 5-NITROTETRAZOLE DETERMINATION OF MERCURY AND COPPER

| Sample | Atomic Absorption | | | | Thiocyanate Titration | | | |
|------------|-------------------|----------------|----------------|----------------|-----------------------|----------------|---------|----------------|
| | Percent Hg(II)(a) | | Percent Cu(II) | | Percent Hg(II)(b) | | | |
| | Found | Ave. of calcd. | Found | Ave. of calcd. | Found | Ave. of calcd. | Percent | Ave. of calcd. |
| 96-8898-80 | 47.1 | 1.1 | 100.7 | 2.3 | 0.012 | 46.79 | 0.05 | 100.0 |
| Modified | 45.2 | 1.0 | 96.6 | 2.2 | 0.004 | 46.60 | 0.28 | 99.6 |
| WOL(M) | 45.2 | 0.7 | 96.6 | 1.5 | 0.034 | 46.70 | 0.05 | 99.8 |
| WOL(P) | 45.4 | 1.1 | 97.0 | 2.4 | 0.038 | 46.70 | 0.05 | 99.8 |
| Batch 2 | 47.3 | 1.0 | 101.1 | 2.1 | 0.028 | 46.37 | 0.05 | 99.1 |
| " 3 | 45.4 | 1.2 | 97.0 | 2.6 | 0.039 | 46.70 | 0.09 | 99.8 |
| " 4 | 45.2 | 0.7 | 96.6 | 1.5 | 0.045 | 46.65 | 0.05 | 99.7 |
| " 5 | 44.8 | 0.7 | 95.7 | 1.6 | 0.037 | 46.79 | 0.09 | 100.0 |
| " 6 | 45.9 | 1.4 | 98.1 | 3.1 | 0.11 | (c) | | |
| " 7 | 46.4 | 1.1 | 99.2 | 2.4 | 0.045 | 45.71 | 0.14 | 97.7 |
| " 9 | 45.8 | 0.7 | 97.9 | 1.5 | 0.042 | 46.46 | 0.19 | 99.3 |
| " 10 | 45.4 | 0.4 | 97.0 | 0.9 | 0.030 | 46.42 | 0.05 | 99.2 |
| " 11 | 45.1 | 0.7 | 96.4 | 1.6 | 0.032 | 47.02 | 0.05 | 100.5 |
| " 12 | 45.3 | 0.8 | 96.8 | 1.8 | 0.046 | 46.51 | 0.09 | 99.4 |

Mercuric oxide plus HNO_3 = reference (see text)


(a) average of two determinations

(b) 96-8898-80, five detns.; Batch 3, three detns.; all others, two detns.

(c) 96.0%, 102.9% (see discussion)

DISCUSSION

For the dichromate oxidation of the 5-nitrotetrazole anion, the equation may be written:

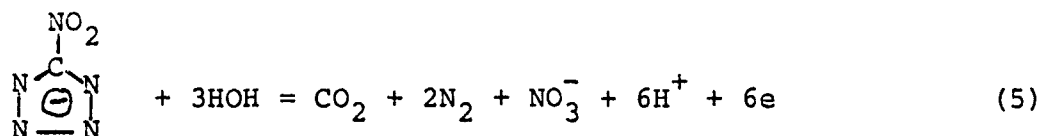

Duplicate analyses for nitrate by a gas chromatographic procedure (6) gave 108% (+ 8%) of the value required by equation (5). Carbon dioxide and nitrogen were not determined. The dichromate oxidation has been used successfully for other explosives in this laboratory (7). The 6 electron change shown above was verified by the experimental results for sample 96-8898-80 of Hg (5-NT)₂, and for Na 5-NT. 2HOH (Table I). The value of 6 was invariant for reflux times of 15 minutes to 4 hours.

Comparing the results for the anion procedures in Table I, the values found by the dichromate procedure are generally 0-3% higher than those found by the UV procedure (the results for Batch #6 will be considered separately below). This may be because the dichromate procedure will measure any readily oxidizable material, for example, organic fibers, or Hg(I).

When the samples of Hg (5-NT)₂ were received, they were under water. A portion was removed and placed in an open container, which was placed in a desiccator over concentrated sulfuric acid to remove any water. To insure that they were dry, the NMR spectrum in deuterated DMSO was observed. No water was found, nor any hydrogen atoms. Thus, organics, such as ethylenediamine, are absent. It is estimated that the NMR method will detect impurities, if present to greater than 3%.

When the results for mercury (II) by the two procedures given in Table II are compared, and the larger average deviation of the AA method is considered, it is seen that the two procedures agree within 0-2%, except for sample 5 where the difference is 2.5%. Generally, the values by AA are lower than those by thiocyanate. When the AA values are higher, it may mean that mercury (I) is present (samples 2, 6, and 7), although there is no experimental proof of this (5).

(6) D. J. Glover and J. C. Hoffsommer, "Gas Chromatographic Analysis of Nitrate and Nitrite Ions in Microgram Quantities by Conversion to Nitrobenzene", J. Chromatography, 94 334 (1974).

(7) For TNT: $\text{C}_7\text{H}_5\text{N}_3\text{O}_6 + 17\text{HOH} = 7\text{CO}_2 + 3\text{NO}_3^- + 39\text{H}^+ + 36\text{e}$; found, 36.1 ± 0.3, RDX: $\text{C}_3\text{H}_6\text{N}_6\text{O}_6 + 9\text{HOH} = 3\text{CO}_2 + 3\text{NO}_3^- + 3\text{NH}_4^+ + 12\text{H}^+ + 12\text{e}$; found, 11.8 ± 0.3

Comparing the data in Tables I and II, the values for the anion by UV appear to be greater than the Hg(II) values by AA, but again, when the deviations are considered, the values by these two methods are in agreement. Except for samples 6 and 7, the anion values by dichromate oxidation agree quite well with the Hg(II) values by thiocyanate. Sample 7 may have Hg(I) present as already mentioned.

Sample 6 is the only sample examined which gives different results for Hg and the 5-NT anion by all four procedures. According to reference (2), during the preparation of sample 6, the desired temperature of 75° could not be reached, and this may account for its low purity.

When treated with ammonium acetate solution, or water only, some of sample 6 did not dissolve. When a large sample of 6 (111 mg) was triturated with a solution of 20% ammonium acetate, diluted with water, and filtered, there was found 7.4% of a yellow solid residue. The filtrate from this solid was analysed for the 5-NT anion by UV and for Hg(II) by thiocyanate titration. 87.7% of the theoretical amount of anion was found and 87.0% of Hg(II). When this percentage of Hg (5-NT)₂, 87.4%, is added to the solid residue, 94.7% of the sample weight is accounted for. This means that 5.3% of the sample weight in solution is not Hg(II) or 5-NT anion.

It is possible that the soluble 5.3% is a Hg(I) compound, because the AA analysis for Hg is 98.1%, and Hg(I) compounds give higher absorptions than Hg(II) compounds (5).


It was also considered likely that a soluble nitrate could be the impurity in sample 6. Consequently, a nitrate analysis was done according to (6). There was found 0.48% nitrate, which corresponds to 0.73% Cu(NO₃)₂. Converting this value of Cu(NO₃)₂ to percent Cu, 0.25%, shows that the nitrate is present not only as the copper salt, for Cu by AA is only 0.11%. Nitrate obviously does not account for the 5.3% soluble in sample 6.

Another possibility for the soluble impurity would be an acidic material. When sample 6 was triturated with water, a pH of 2.99 was found, while the pure sample 96-8898-80 gave a pH of 3.27. When these pH values are converted to H⁺ concentrations, there is found for sample 6, 10.23×10^{-4} M, and for 96-8898-80, 5.37×10^{-4} M. The difference, 4.86×10^{-4} M, would be the normality of the acidic impurity. If the acid impurity is hydroxytetrazole (molecular weight = 86), then the weight of hydroxytetrazole is $86 \times 4.86 \times 10^{-4} = 0.0418$ g/liter. As the sample weight was 0.006467 g/10 ml, hydroxytetrazole content would be 6.5%. This value is in fair agreement with the experimental value of 5.3%.

Titration of sample 6 with thiocyanate without filtration of the insoluble material, gave a first end point at 96% of the calculated Hg(II). This end point faded, and a final end point was found at 102.9% of the calculated value. Since the end point did not fade

when the insoluble material was filtered, the insoluble material is probably a Hg(II) compound. Considering the value of 102.9%, the molecular weight of the insoluble material must be less than that for Hg (5-NT)₂, 428.7 g/mole. It is noted that mercuric 5-hydroxy-tetrazole has a molecular weight of 370.6 g/mole.

The dichromate oxidation of sample 6 without filtration of the insoluble material, gave a value of 94.2% of the sample weight, based on the expected 6 electron change for 5-nitrotetrazole. If, indeed, hydroxytetrazole is present, either as the acid or Hg(II) salt, then the electron change expected is 4:

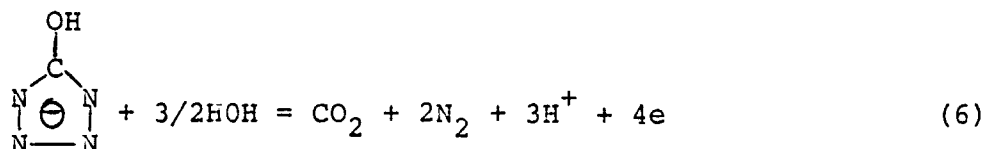

This means hydroxytetrazole as an impurity would give low results by the dichromate oxidation. However, these results will be higher than the anion by UV, if the 5-hydroxytetrazole anion absorbs less than the 5-NT anion, which is what is observed (Table I).

SUMMARY AND CONCLUSIONS

Mercury (II) in mercuric 5-nitrotetrazole may be determined by atomic absorption using the unique-to-mercury absorption line at 2536.5 angstroms, with a precision of 1-3%. For a precision of 0.1 to 0.6%, a titration with thiocyanate ion is recommended. The 5-nitrotetrazole anion may be determined by measuring the absorbance due to the anion at 257 nm in the ultraviolet, with a precision of $\pm 1\%$. The anion may also be determined by oxidation with dichromate ion in sulfuric acid solution, with a precision of $\pm 0.5\%$.

For specification purposes, there are advantages to using both of the cation methods presented here, as well as using both anion methods. The AA method for Hg depends on measuring the unique-to-mercury absorption line at 2536.5 angstroms. Hg(I), if present either from Hg(II) contamination, or by reduction of Hg(II), will show up by giving higher than calculated values based on Hg (5-NT)₂. The thiocyanate titration gives more precise values for Hg(II), and will react with Hg(I) to give the insoluble salt, and thereby give low results. The UV method for the anion is specific, but is subject to interference by any material absorbing at 257 nm. The dichromate oxidation of the anion will also reflect any easily oxidizable impurity.

Examination of the results for the fourteen different samples of Hg (5-NT)₂ presented in Tables I and II, show that one (Batch 6) gives analyses indicating a poor product. On the basis of all the data, the following requirements are recommended:

| <u>Property</u> | <u>Requirement</u> |
|--|--------------------|
| Assay | |
| Hg (AA) | 100 \pm 3% |
| Hg (thiocyanate) | 99.5 \pm 0.5% |
| 5-nitrotetrazole anion | |
| UV | 100.0 \pm 1.0% |
| Dichromate oxidation | 100.0 \pm 1.0% |
| Cu (AA) | <0.05 \pm 0.01% |
| NO ₃ ⁻ (ref. 6) | <0.10 \pm 0.01% |
| Insoluble in 20% ammonium acetate | 0.5% max. |
| pH (10 mg triturated with 10 ml water) | 3.32 \pm 0.10 |

DISTRIBUTION LIST

Copies

| | |
|-------------------------------------|---|
| Commander | |
| Naval Sea Systems Command | |
| Washington, D. C. 20362 | |
| SEA-0332, A. Amster | 1 |
| SEA-0332B, G. Edwards | 1 |
| SEA-04H, E. Daugherty | 1 |
| SEA-04JE, R. Beauregard | 1 |
| SEA-09G32, Technical Library Branch | 2 |
| Commander | |
| Naval Air Systems Command | |
| Washington, D. C. 20361 | |
| ESA 2041 D1, H. Graham | 1 |
| Chief of Naval Material | |
| Washington, D. C. 20360 | |
| Code 0323, I. Jaffe | 1 |
| Commanding Officer | |
| Naval Weapons Support Center | |
| Crane, Indiana 47522 | 1 |
| Commanding Officer | |
| Naval Ammunition Depot | |
| c/o Fleet Post Office | |
| San Francisco, California 96612 | |
| Library | 1 |
| Commanding Officer | |
| Naval Ordnance Station | |
| Indian Head, Maryland 20640 | |
| Library | 1 |
| Superintendent | |
| Naval Post Graduate School | |
| Monterey, California 93940 | |
| Director | |
| Naval Research Laboratory | |
| Washington, D. C. 20390 | |
| Technical Information Section | 1 |
| Commander | |
| Naval Weapons Center | |
| China Lake, California 93555 | |
| Technical Library | 1 |
| E. Royce | 1 |
| I. Joyner | 1 |
| B. Hayes | 1 |

DISTRIBUTION LIST (Continued)

Copies

Commanding Officer
 Naval Weapons Station
 Yorktown, Virginia 23491
 Library
 W. McBride
 L. Rothstein

1
 1
 1

Office of Naval Research
 Department of the Navy
 Arlington, Virginia 22217
 R. Miller

1

Naval Medical Research & Development Command
 National Naval Medical Center
 Bethesda, Maryland 20014
 LCDR Leigh Doptis, Code 47

1

Naval Medical Research Institute
 National Naval Medical Center
 Environmental Biosciences Department
 Bethesda, Maryland 20014
 LCDR Thomas A. Hill

1

Commanding General
 Air Force Systems Command
 Andrews Air Force Base
 Camp Spring, Maryland 20334
 Library

1

Commanding General
 Army Material Command
 Department of the Army
 Washington, D. C. 20318
 R & D Division

1

Acting Director
 Army Material Systems Analysis Agency
 Building 392
 Aberdeen Proving Ground, Maryland 21105
 Library

1

DISTRIBUTION LIST (Continued)

Copies

| | |
|--|--------|
| Commanding General Army Missile Command Redstone Scientific Information Center Redstone Arsenal, Alabama 35809 Technical Library | 1 |
| Commanding General Aberdeen Proving Ground Aberdeen, Maryland 21001 Technical Library | 1 |
| Commanding Officer U.S. Army Armament Research and Development Command Dover, New Jersey 07801 Library | 1 |
| R. Walker | 1 |
| H. Matsuguma | 1 |
| J. Hershkowitz | 1 |
| W. Voreck | 1 |
| Commanding Officer Harry Diamond Laboratories 2800 Powder Mill Road Adelphi, Maryland 20783 DRXDO-TI, Scientific & Technical Information Office K. Warner | 1 1 |
| Commanding Officer Eglin Air Force Base, Florida 32542 ADTC (DLOSL), Technical Library | 1 |
| Commander Air Force Armament Laboratory RTD Weapons Division Technology Branch (ADLRW) Eglin Air Force Base, Florida 32542 Library | 1 |
| Defense Documentation Center Cameron Station Alexandria, Virginia 22314 TIPCR | 12 |

DISTRIBUTION LIST (Continued)

Copies

Director
 U.S. Bureau of Mines
 Division of Explosive Technology
 4800 Forbes Avenue
 Pittsburgh, Pennsylvania 15213

1

Director
 Los Alamos Scientific Laboratory
 P.O. Box 1663
 Los Alamos, New Mexico 87544
 Library GMS-7

1

NASA
 Langley Research Center
 Hampton, Virginia 23365
 Library

1

NASA
 Lyndon B. Johnson Space Center
 Houston, Texas 77001
 Library
 T. Graves
 M. Falbo

1

1

1

NASA Scientific & Technical Information Facility
 P.O. Box 33
 College Park, Maryland 20740

1

Lawrence Livermore Laboratory
 University of California
 P.O. Box 808
 Livermore, California 94551
 Technical Information Division
 M. Finger

1

1

Allegany Ballistics Laboratory
 Cumberland, Maryland 21501
 Library

1

Jet Propulsion Laboratory
 California Institute of Technology
 Pasadena, California 91125
 Library

1

DISTRIBUTION LIST (Continued)

| | Copies |
|--|--------|
| Stanford Research Institute Propulsion Sciences Division Menlo Park, California 94025 Marion Hill | 1 |
| R. H. Stresau Laboratory Star Route Spooner, Wisconsin 54801 | 1 |
| Thiokol Chemical Corporation Redstone Arsenal Huntsville, Alabama 35809 | 1 |